
Onyx® Vantage 9590



"The Onyx unit I use has been dropped on hard floors, hard rocks and has been splashed with cold water. The Onyx continues to work in all these environments and remains very readable and usable."

- Carl R., RN

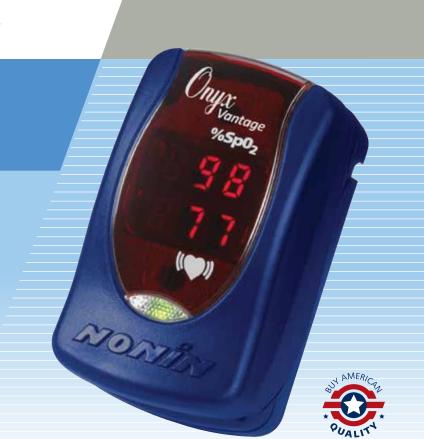

Accuracy you can act on™.

When you need to make decisions quickly, you do not have time to question the accuracy of your finger pulse oximeter. You can trust Nonin Medical's $Onyx^{\circ}$ Vantage finger pulse oximeter with $PureSAT^{\circ}$ technology to quickly and accurately capture SpO_2 and pulse rate measurements — even on patients where low perfusion is a challenge. In fact, Onyx is the only finger pulse oximeter brand whose accuracy and performance claims are supported by published peer-reviewed, clinical studies. Accuracy matters, and Onyx Vantage provides accurate data for actionable decision making.




Onyx® Vantage 9590

Product Highlights

The Onyx is the only finger pulse oximeter with scientifically proven accuracy and is tested for use on fingers, thumbs and toes.

- **Proven Accuracy** Scientifically proven accuracy in the widest range of patients.
- Reliable Results Accurate readings on your challenging cases due to low perfusion and darker skin tones. And, Onyx Vantage 9590 has been
- Versatile One unit works on pediatric to larger adult patients and can be used on fingers, thumbs or toes.
- Cost Efficient 4-year warranty and up to 6,000 spot checks on

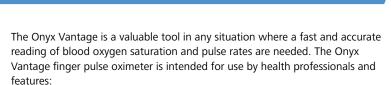

Clinician-facing, easy-to-read, bright LED – SpO₂ and pulse rate can be read from any angle, in all lighting conditions, day or night.

Three-color pulse quality indicator LED - Green, yellow or red light provides a quick assessment of patient pulse quality.

Accuracy in the widest range of patients and conditions – Only Nonin uses PureSAT® SpO₂ technology – an intelligent pulse-by-pulse filtering to provide accuracy in situations that are challenging due to dark skin tones, low perfusion and more.

One product for pediatric to adult patients -

Automatically adjusts to each patient from pediatric to adult quickly and accurately and accommodates fingers, thumbs and toes with digit thicknesses: 8 mm to 25.4 mm (0.3" to 1.0").





Proven durability – Onyx Vantage is tough enough to withstand a minimum of 50 drops, is rugged enough to be used in the most demanding situations and exceeds IP32 water ingress testing.

Automatic turn on and off - No on/off button! Just insert digit for fast and easy spot checks. The auto-on/off mechanism allows two AAA batteries to provide up to 6,000 spot checks before being replaced!

Portable – Take the Onyx Vantage anywhere — it is lightweight and weighs less than two ounces.

Available in four colors - Black, Red, Blue and Purple. Complete with lanyard, instructions on CD and two AAA batteries. Optional carrying case and retractable lanyard available.







"I've learned that the accuracy and reliability is quite good. I've managed to use it on patients that perhaps were not getting a peripheral saturation using other sources or saturation monitors, but I usually could get it on my little Onyx fingertip from Nonin." M. Shaw, RT

Specifications

Dimensions	5.59 cm (2.2") L x 3.3 cm (1.3") W x 3.23 cm (1.27") H
Warranty	. 4 years
Oxygen Saturation Display Range	0 – 100% SpO ₂
Pulse Rate Display Range	18 – 321 beats per minute (BPM)
Oxygen Saturation Declared Accuracy Range (A _{rms} *)	70 – 100% SpO ₂ ±2 digits
Low Perfusion Oxygen Saturation Declare Accuracy Range (A _{rms} *)	
Pulse Rate Declared Accuracy Range (A _{rms} *)	20 – 250 BPM ±3 digits
Low Perfusion Pulse Rate Declared Accura Range (A _{rms} *)	
	ower** . 660 nanometers @ 0.8 mw maximum averag . 910 nanometers @ 1.2 mw maximum averag
Temperature Operating	
Humidity Operating Storage/Transportation	
Altitude	-
Operating Hyperbaric Pressure	
Battery Life Operating	. Approximately 6,000 spot-checks or 36 hours of continuous operation using new AAA-size
Storage	alkaline batteries . 12 months
Classifications per IEC 60601-1 / CAN/CSA-C2: Degree of Protection	
Mode of Operation	. Continuous

This product complies with ISO 10993-1, Biological evaluation of medical devices -Part 1: Evaluation and testing.

- *±1 A_{rms} represents approximately 68% of measurements.
- **This information is especially useful for clinicians performing photodynamic therapy.

Nonin Medical B.V.

Prins Hendriklaan 26

Tel: +31 (0)13-79 99 040

Fax: +31 (0)13-79 99 042

E-mail: infointl@nonin.com

1075 BD Amsterdam • Netherlands

Specifications are subject to change without notice.

Nonin Medical, Inc.

13700 1st Avenue North

Plymouth, MN • 55441-5443 • U.S.A.

Tel: +1.763.553.9968 1.800.356.8874

Fax: +1 763 577 5521

E-mail: info@nonin.com

onyxvantage.com





